


CLEANTYPE[®] PRIME SERIES

PANEL+ (WITH BACKLIGHT)
PANEL

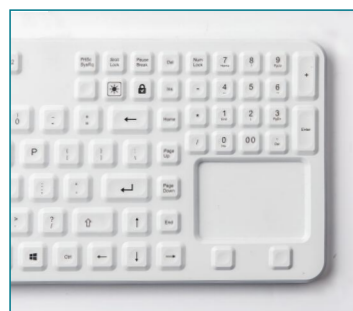

SPECIAL FEATURES

- Adjustable brightness of backlight (PANEL+)
- 109 keys incl. num pad and track pad
- Match with surgical panel (wall mount)
- Excellent tactile feeling
- Integrated Clean Mode for cleaning
- EN60601-1-2 Ed 4 certified

SPECIFICATION

Number of keys:	109 (US-Layout)	Interfaces:	USB 2.0
Keyswitch technology:	carbon contact	Cable length:	15cm
Switching cycles:	approx. 2 mill. (per key)	Operating temperature:	0°C ~ +55°C
Housing material:	silicone	Storage temperature:	-30°C ~ +70°C
Housing color:	white & black	Dimensions:	365(L) x 134(W) x 14.2(T) mm
Backlight:	blue color, 5 level adjustable brightness	Weight:	approx. 415g
		Certification:	EN60601-1-2 Ed 4, CE, FCC, RoHS, REACH
		Warranty:	2 years


ORDER INFORMATION

Cat. No.	Product description	Color	Layout
KSI-U10220	CLEANTYPE [®] PRIME PANEL+	White	US/DE/JP/KR/THAI
KSI-U10221	CLEANTYPE [®] PRIME PANEL	White	US/DE/JP/KR/THAI
KSI-U20220	CLEANTYPE [®] PRIME PANEL+	Black	US/DE/JP/KR/THAI
KSI-U20221	CLEANTYPE [®] PRIME PANEL	Black	US/DE/JP/KR/THAI

Other layouts, configurations and interfaces on request.

Customized key function on request (Firmware update via USB)

TECHNICAL DRAWING (IN MM)

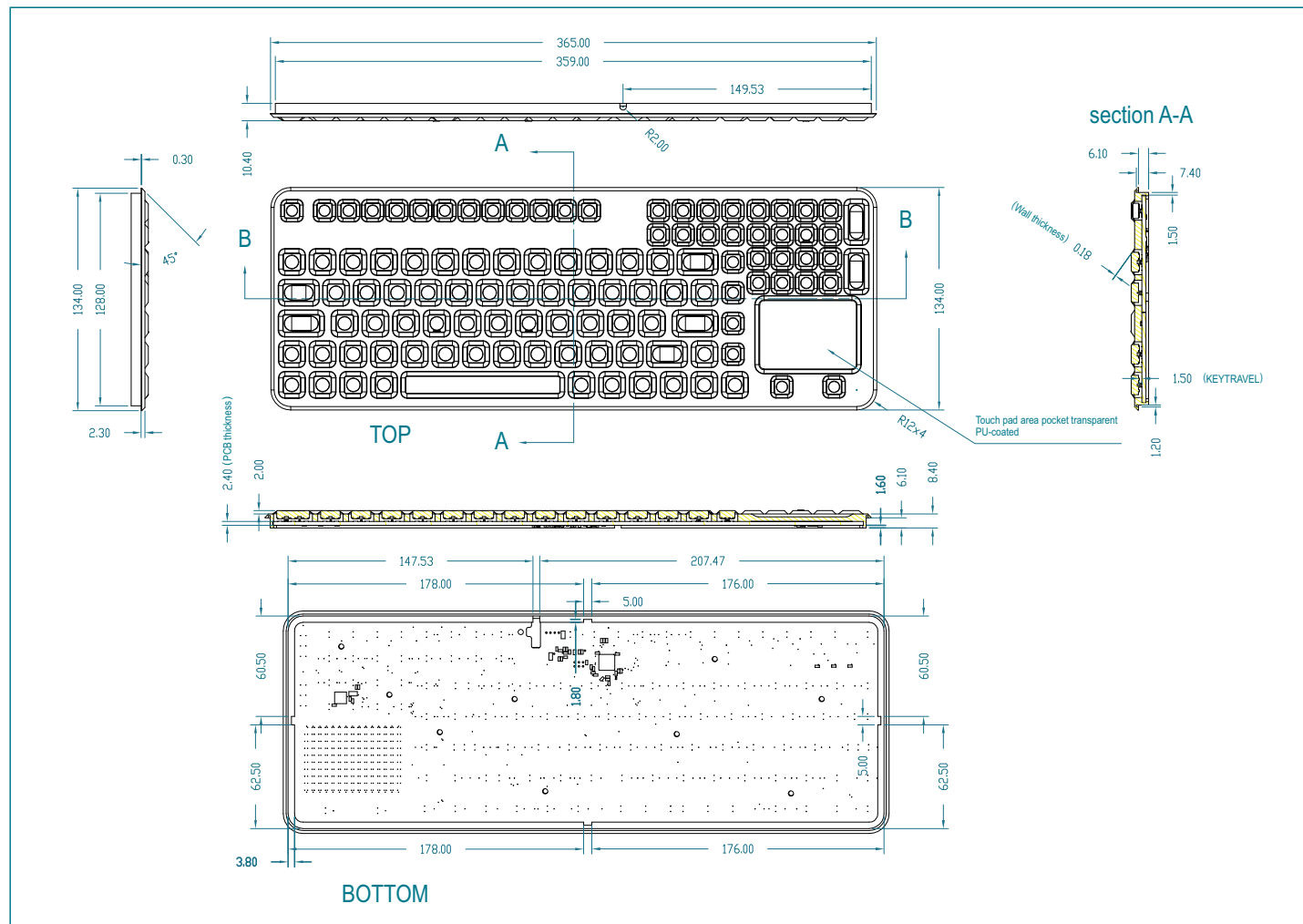

EN60601-1-2 Ed 4 certified

EN 60601-1-2 Ed 4 applies to the basic safety and essential performance of Medical Equipment(ME) equipment and ME systems in the presence of electromagnetic disturbances and to electromagnetic disturbances emitted by ME equipment and ME systems.

SHORT OPERATING INSTRUCTIONS



Press and hold lock key for 5s to enable Clean Mode.
Clean Mode is indicated by breathing backlight.



Adjust backlight level with this button.
Five level can be adjusted including off.

CLEANING INSTRUCTION

For further details please refer to below link:
<http://www.gett-group.com/content/service-cleaning-instructions>